



EXPERIENCE THE DIFFERENCE
YOUR SECURITY & SAFETY PRODUCTS

RAPID TEST For Coronavirus Disease



RAPID TEST Antibody (IgG/IgM) COMBINED TEST KIT



- RAPID Results within 5-10 minutes
- DIVERSE Serum, Plasma, Whole blood
- EASY Only 4 steps taken
- RELIABLE IgG/IgM tested at one time




Specification:


Product name	Coronavirus Disease Antibody (IgM/IgG) Combined Test Kit
Test principle	Lateral flow immunoassay
Sample	Whole blood, Serum, Plasma
Time to report	5-10 minutes
Operation temperature	Room temperature
Pack	25 tests
Material provided	25 individual sealed pouches, 25 disposable droppers

RAPID TEST IgM/IgG antibody combined rapid test for COVID-19 is used to qualitatively detect IgG and IgM antibodies of the novel coronavirus (SARS-CoV-2) in human serum, plasma or whole blood in vitro. These biomarkers indicate that a person's immune system has responded to infection with the COVID-19 virus.

This innovative rapid test kit is developed with lateral flow immunoassays technology in high sensitivities and specificities. It is quite easy to use and ideal for screening check in epidemics control of COVID-19 virus.

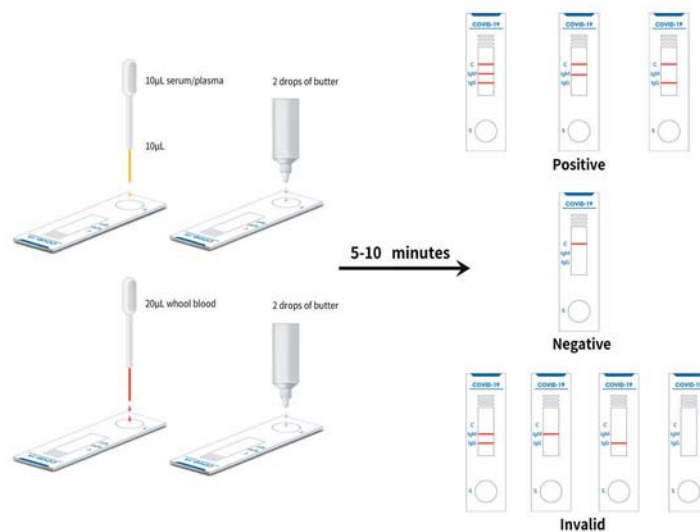

- 1 Collect blood sample
- 2 Place 2 drops of buffer in sample well
- 3 Add 10 pL serum/plasma or 20 pL whole blood to sample well
- 4 Read results in 5-10 minutes

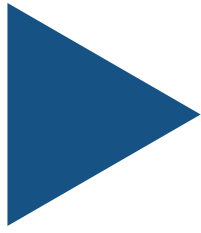

LEARN MORE: c4icommunication.com

FOR MORE INFORMATION: contact@c4icommunication.com or +1 (302)9811340